

## Reliable Quality

The ePM 12M Vet comes equipped with a high-quality shell and accessories. They are safe, durable, and built to last.

### High-quality shell materials

- Manufactured with industry-leading plastics that stand up to the most stringent medical-grade cleaners.

### Fan-free design

- Keeps fur out of equipment
- Reduces cross-contamination, protecting the health of patients and staff
- No internal dust, extending the life of the monitor

### Veterinary-specific accessories

- Low-trauma ECG clips maintain consistent contact
- Appropriately-sized NIBP cuffs provide accurate readings
- Durability-tested for long-lasting reliability



ECG: Non-traumatic alligator clips



NIBP: Variety of cuffs fit multiple patient sizes

## ePM 12M Vet

A new generation of veterinary modular monitor

## Intelligence for safety



As veterinary medicine advances, the demand for high-quality veterinary equipment increases. The Mindray ePM 12M Vet focuses on progressive design, making this monitor smart and simple.



- High-Definition Capacitive Touch Screen
- Gesture controls creates a smartphone-like user experience
- Two-step operations streamline use
- Wireless capability

## Excellent Performance

ePM 12M Vet supports advanced parameters with excellent and accurate performance.

### Precise algorithms

- NIBP: Veterinary-specific; weak-signal tolerant and anti-interference algorithm
- ECG: Patented Mindray exclusive multi-lead analysis algorithm
- SpO<sub>2</sub>: Masimo technology

### Comprehensive parameter monitoring

- Innovative design of lifting handle and flexible expansion
- Supports advanced parameters: CO<sub>2</sub>, 4-IBP, C.O., AG, O<sub>2</sub>
- CO<sub>2</sub> and A.G. : Artema technology, with 30 years of technical expertise.
- Up to 3-channel recording and printing




## Smart & Simple

The ePM 12M Vet provides doctors and nurses with smart applications and richer intraoperative information, making monitoring more convenient, and supporting surgery more efficiently.


### Smart lead-switching & Visible alarm

- Automatically switches to an available lead if the main ECG lead detaches, ensuring continuous monitoring.
- One click confirms which lead is detached, increasing troubleshooting efficiency

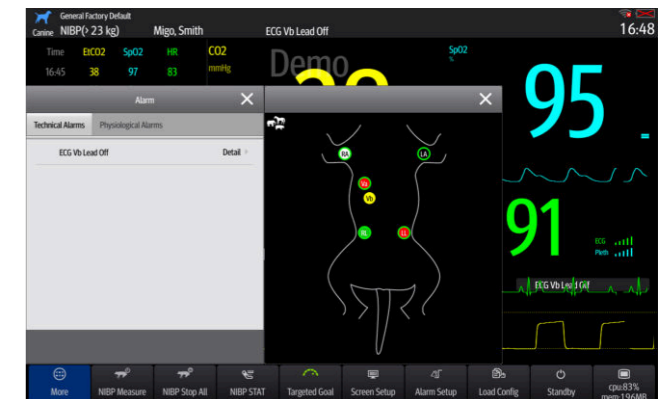

### Target reference range display

- Target range dashboard for key parameters helps focus on key parameters during surgery
- Customizable target ranges and customizable key parameter selection

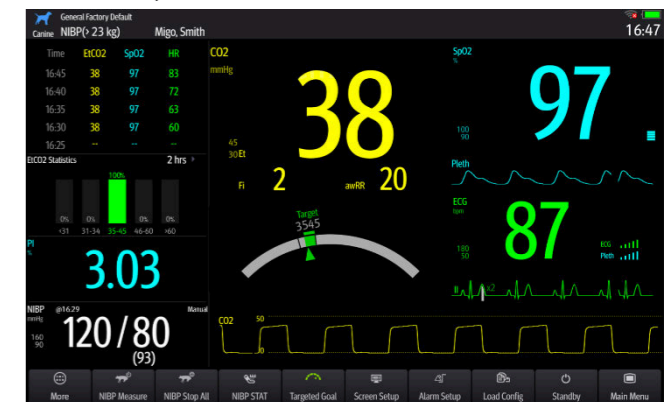

### Customized configuration

- Memory for up to 25 customized configurations. Identifiable configurations with different alarm limit settings for different species or patients
- Supports one-click USB import/export – configurations can be loaded to other ePM monitors

