

Wireless Color Doppler Ultrasound Scanner

For bovine, equine, ovine scanning

S10

BestScan®



Professional Veterinary Applications

BestScan® S10

Greater accuracy & earlier detection compared. Suitable for Bovine, Equine, Ovine, Swine, widely used in uterus, ovaries, fetal structure, fetal sex, reproductive tract disease, horse's far limb, tendon, bone, eyes, abdomen and chest organs.



Platform Features

High performance

AI intelligence

High speed

Stability

Security

Application platform

Android system, Window system, ISO system
Support platform

Viewing device

Smart phone; Tablet PC; VR Goggles

Android




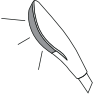

Android 4.4 (API 19) or above, equipped with WI-FI
802.11n and BLE 4.1 models * support blue tooth.




Multi-language interface

supports 8 languages:
Chinese-English-Arabic-German-Spanish-French-Italian-Russian

Unique Design for Veterinary

- ▶ High-end atmosphere, sealed dust and waterproof integrated design, IP67 waterproof performance, From the probe tip to the display, easy to clean.
- ▶ Only 650 g with compact, lightweight design, easy to carry.
- ▶ Double-sided anti-slip for easy grip. ABS+TPU material, strong wear resistance and crash-proof. (drop test of 1 meter) the probe cable is flexible and bendable, suitable for the difficult environment of the farm.
- ▶ Button finger touch recognition for blind operation, Separate freeze button design . 6 function keys adjust the gain, depth, and freeze separately. Sensitive keys and easy to operate.

 <p>HD image</p>	<p>Broad-bandwidth patented probe rectal linear probe and rectal convex, getting images from deeper and wider direction. Excellent image quality helps identification of early pregnancy, fetal gender and uterine diseases easily and early.</p>
 <p>Accurately</p>	<p>Accurately measure blood flow velocity, detecting the hemodynamic status of blood vessels, increasing Independently designed & patented for probes, high quality image performance .</p>
 <p>Penetration</p>	<p>Probe penetrating enhanced, getting images from deeper and wider direction.</p>
 <p>Image field</p>	<p>High sensitivity sensor for better images, unique design to improve the accuracy of late pregnancy diagnosis, especially C-type probe is very light, and the probe has a large open space on one side and a large operable area.</p>
 <p>Frame rate</p>	<p>Image processing with AI, rapid imaging and high quality, 26 frame rate / second image, when fast moving probe, no delay or lag.</p>

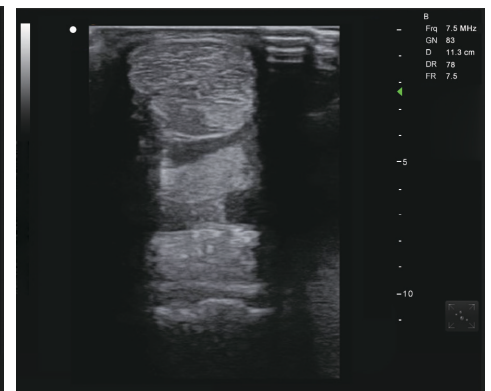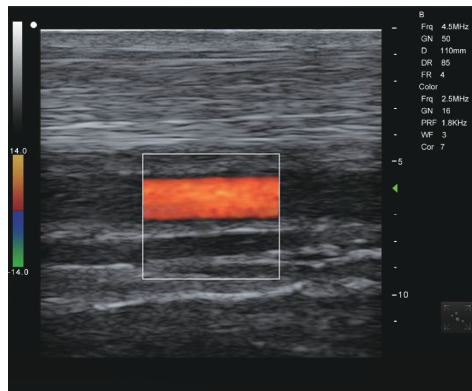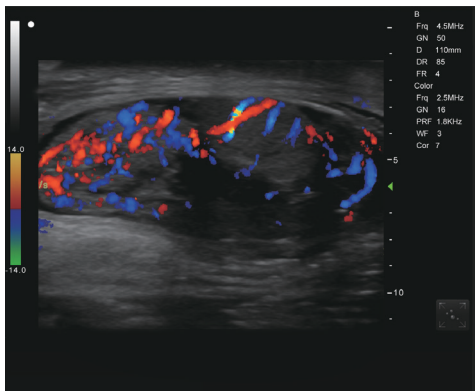

BestScan® S10



The S10 app is available on compatible Android smart devices via scanning code to download



Fast Connection and Conversion

BestScan® S10



Supports connectivity to multiple smart devices such as PAD, Mobile phone, laptop (Android, IOS, Windows) via WI-FI in 10 S within fast connection; Tablets and TVs can be connected at the same time, Supports USB-DICOM, Wireless printers, HDMI.

Power to go, wireless charging, battery built-in integration 12 hours working time, video glasses support external mobile power for up to 4 hours.



Rich measurement software packages: professional vet measurement, Grids for management, measurement software packages to calculate gestational age, expected date of birth, etc.

Simple and Fast Work Workflow

Easy to Use

Classic 3-step
workflow

No connection
password

Second intelligent
identification password



Step1: Prepare machine



Step2: Farm scanning











Step3: Cleaning

Starting up -Open intelligent terminal software - Connect - Start diagnosis. Through the tablet and cable box converters to achieve wired multi-screen interaction (up to three screens)

Plenty of Accessories Options

BestScan® S10

	<p>I-scan smart goggle Customized Binocular Type HD Wearable Video Glasses</p>
	<p>Rectal Probe Introducer Rounded 85cm,eases your shoulder and arm strain for Fetal Aging, Pregnancy Testing</p>
	<p>Chest Mount for Tablet To be used in conjunction with your choice of tablet</p>
	<p>Wrist Mount for phone To be used in conjunction with your choice of Phone</p>

	<p>Field working bag Storage the probe and its accessories</p>
	<p>Ipad with S10 software installed Use as a viewing device (manufactured by Apple Inc)</p>
	<p>Power bank 10000mAh, 147x71.2x14.2mm, supply power for goggle</p>
	<p>Mobile Car Charger Battery Charging in the car</p>

After Sales Service

Wireless upgrades

Full 2-Year Warranty for main unit

An affordable teaching tool
Professional veterinary technical support

ADD: 4/F, BMV of Yinjin Building, Block 71, Baoan Centre District, Shenzhen 518010, P.R. China

Tel: +86-755-26564580

E-mail: sales@bmv.cc

Web: www.bmv.cc

BMV is first brand in the veterinary industry of China.

Set R&D, production, sales in one of the medical equipment manufacturers, service providers

BMV & BestScan® are registered trademarks or trademarks owned by BMV Technology Co., Ltd.

©2008-2019 BMV Technology Co., Ltd. All rights reserved.

Specifications subject to changes without prior notice.

P/N:ENG-S10-BROCHURE-210285X1 0P-201906

